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G4+

PLABEAU G4+ User Manual

* Ensure you familiarize yourself with the user manual prior to use.

PLABEAU-Certified by Korea
Asthma Allergy Foundation
and Korea Atopy Association.
CE & CCC labelling existing.

CONTENTS

01	Compositional items and names of parts	20
02	Operating method / How to use	22
03	Charging method / Explanation on LED	24
04	FAQ / Product specification	25
05	Precautions for users	26
06	Warranty Terms & Conditions / Others	28


Plasma beauty solutions

PLA BEAU G4+


Balance and restore your skin with just 10 minutes of PLABEAU a day.

PLABEAU utilises innovative cold plasma technology to engineer state of the art solutions in multiple areas of skin management. From anti-aging and deep cleansing to overall skin balancing, PLABEAU is meticulously designed to provide your ultimate skin care experience.

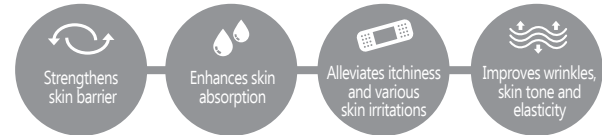
01 Compositional items and names of parts

Components: Images of the components displayed below may differ from the actual components and are subject to change in accordance to changing circumstances.

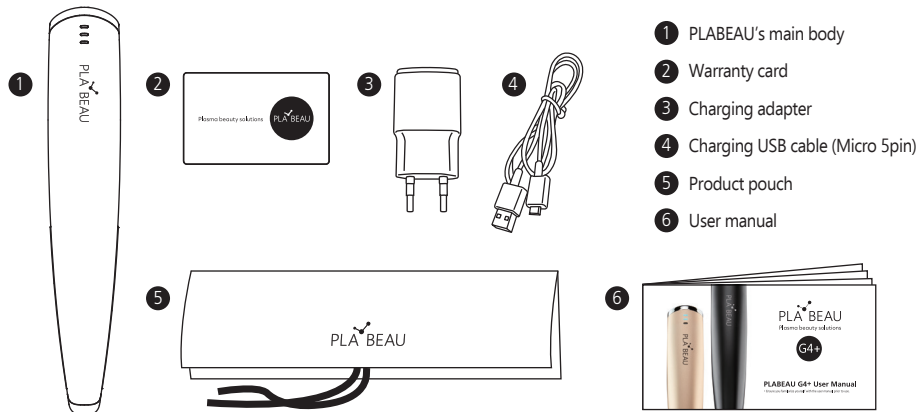Name of parts and components :

Power and LED mode

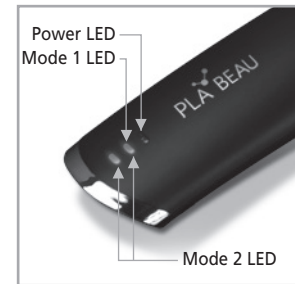Ceramic surface pad
emits plasma ions

Power and mode switch
Mode 1: Sensitive Mode
Mode 2: Normal Mode

Charging cable connection part
(Please use a recommended charger.)

Operating method of PLABEAU

STEP 1

Please turn the device ON/OFF by sliding the power button located on the center of the device. When the power button is on Mode (2), the device operates on Normal Mode, and when the power button is on Mode (1), the device operates in Sensitive Mode.

! When PLABEAU is on for 10 minutes, a blue LED light flickers, indicating that the device has stopped generating plasma. Flickering blue LED light indicates recommended time of use.

STEP 2

Please connect the charging port to the main body when Red LED light blinks. Red LED light is on when it is charging.

! Within PLABEAU's built-in battery, user can use PLABEAU on broad skin parts 4 ~ 5 times. (app. 1 hour)

STEP 3

Green LED light turns on when battery is fully charged.

How to use PLABEAU 1


Please completely wash your face and other areas of the skin you intend to apply PLABEAU on. If the residue of make-up and other cosmetic products are left on the skin, it may reduce the efficacy of the device. Please make sure to use PLABEAU when your skin is completely dried.

! If water is left on the skin, it may cause a stinging sensation. Please use this product after your skin has completely dried.

How to use PLABEAU 2

Please scan and massage the device on your skin while holding onto the grip section of the device. If you want more intensive care, set the device on Mode (2): Normal Mode and apply PLABEAU on areas you intend to focus on for 5 to 10 seconds.

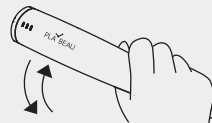
! We recommend applying the device for 10 minutes a day. Depending on your skin sensitivity, the level of stinging you experience may be different. If you have sensitive skin, start by using the device for 2 to 3 minutes and incrementally increase usage of the device.

! When plasma is generated, a blue LED light will turn on. It is normal to hear noises and smell a metallic odor as these are all byproducts of generating plasma.

How to use PLABEAU 3

Please apply moisturizers on your skin after the use of PLABEAU. Before applying the moisturizers, you can clear up your skin texture by lightly washing off face or applying toner over your face.

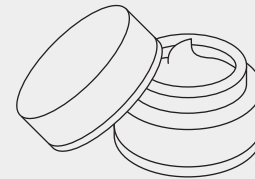
How to use PLABEAU 4

After the use of PLABEAU, sebum and skin impurities may remain on the plasma ceramic pad. We recommend keeping the plasma pad clean by regularly swiping it with alcohol and a cotton pad after use.


03 Charging method / Explanation on LED


1 Components for the charging process.

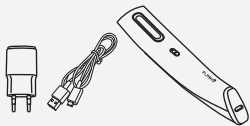
3 Insert the plug of the charging adaptor into the power outlet.

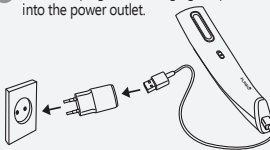2 Connect the charging jack with the main body.

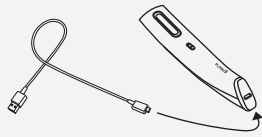4 Once charging is completed, separate the product and adaptor before pulling out the plug from the power outlet.

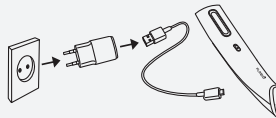LED Indicator and Operating Time

• Power on:

- ① Sensitive Mode (Switch 1): Middle blue LED light turns on
- ② Normal Mode (Switch 2): Top and middle blue LED lights turn on

- 10 minutes after operation: Blue LED blinks
- During recharging: Red LED
- Fully Charged: Green LED Light
- When recharging is needed: Red LED blinks

- * When Red LED flickers, recharge the device immediately.

- Device will not operate during charging.
- Recharging the unit with unapproved recharger will cause explosion or breakdown of the product.
- It is fine to use the commercially available Micro USB (5-pin) connector type recharger for recharging.
- Product may breakdown if the recharger is connected erroneously.
- Detach recharger when not using the unit to conserve energy.

- Product warrantee will not cover the damages arising from erroneous use.
- You must always disconnect the power outlet and recharger to prevent wasting of electric power.
- Make sure to fully recharge the unit when you are using the product for the first time after the purchase or after having left the product unattended period of time.
- Keep out of reach of children.
- Device is not waterproof!

04 FAQ / Product specification

Q I want to improve my skin appearance as quickly as possible. How often should I use the PLABEAU device?

A Although it is recommended that you use PLABEAU twice a day, 10 minutes per session. When you start using the device, you should use it for a shorter period of time (2~3 minutes) and gradually increase usage to allow time for your skin to adapt to the device.

Q I smell metallic odor while using PLABEAU. Is it normal?

A Smelling a metallic odor is completely normal and is a natural byproduct of generating plasma.

Q I experience a warm sensation on my skin when using PLABEAU.

A It is normal to experience warmth while using PLABEAU. Difference in sensation may vary depending on the sensitivity of your skin. If you experience any form of inconvenience due to the warm sensation, stop use and consult with our Customer Center.

Q Skin becomes reddened or trouble occurs after the use of PLABEAU.

A We recommend that you stop using the device for about 2 to 3 days before recommencing. If reddening or irritation does not reduce, consult a dermatologist before using the device again.

Product name	PLABEAU G4+
Product code	PB-G4P
Area of care	Entire face and body
Effectiveness and efficacy	Cosmetic effects on skin and promotion of cell regeneration
Power source	DC 5V, 1.5A
Electric power for plasma	less than 0.75W
Material	Aluminum & Polycarbonate (PC)
Size	35(W) x 200(L) x 37(H)mm
Weight	Approximately 120g
Certification No.	MSIP-REM-PBZ-PB-G4P
Manufacturer	Plabio Co., Ltd.
Made in	Korea

Use only as intended!

05 Precautions for users

Before using this product, be sure you are familiar with the precautions.

-  **Precaution 1** Patients exhibiting the following symptoms should use the product only after consultation with a medical specialist:
1. Those with skin rash or experiencing irritated or stinging skin.
 2. Those with inflammation of the skin.
 3. Those with particularly sensitive skin.
 4. Pregnant women.
 5. Those taking medicines that may cause skin sensitivity.

-  **Precaution 2**
1. Do not use the product for purposes other than its recommended purpose.
 2. Do not use products by inserting numerous electrical products into the same power outlet.

-  **Precaution 3** Users of the following medical devices are not allowed to use PLABEAU:
1. Users of electronic medical devices that are worn inside the body such as cardiac pacemaker.
 2. Users of electronic medical devices for sustaining of life such as artificial cardiopulmonary machine, etc.
 3. Users of electronic medical devices attached to the body such as electrocardiogram (ECG) device, etc.

*Note: Those with sensitive skin may experience slight skin trouble upon contact with the device.
If symptoms persist, cease using the device for 2~3 days before recommencing.

 Highly dangerous situations

Do not put in the mouth

Do not let the device come in contact with eyes

Do not use with wet hand

Do not use the unit if the section for generation of plasma is cracked or broken

Never use this device in the bathtub, shower or sink filled with water

Please ensure to abide by the precautions when using the device

1. Never use directly on eyes.
2. Make sure to take off all jewelry from your body before using the device.
3. Make sure during use that the plasma pad does not crack or break, otherwise there is a risk of skin injuries, as burns can occur. Compensation for negligently caused damage is not possible.
4. Never store or use in places with high temperature and humidity, such as saunas or bathrooms.
5. Never pull the cable out with a wet hand as you risk electrocution.
6. Familiarize yourself with the instructions for use and use the device only in a dry location and use the device appropriately.
7. Stop using the device if you detect any abnormalities during use and turn off the power device off. Contact your dealer immediately.
8. Avoid situations where the device may malfunction, e.g. when in contact with other electrical equipment.
9. In the event of thunderstorms, interrupt the application and remove the device from the plug.

Please ensure to abide by the precautions in storage

1. Do not store the product in locations with high humidity or temperature and in very dusty places with direct exposure sunlight. Store in a dry place!
2. Avoid storage in locations with magnetic or electromagnetic waves.
3. Keep the device out of the reach of infants, children, people with dementia and pets.
4. Do not drop the device in the water. A free repair is not possible for negligent damage.
5. Always use and store the device under hygienic conditions. Foreign matter on the product that will not be removed can lead to a failure.
6. Never store the plasma emitting ceramic surface (plasma pad) on top of it.

Precautions regarding maintenance

1. Do not disassemble the device and do not attempt to repair it yourself. In case of failure, interrupt they use and contact the dealer.
2. Do not drop the unit or subject it to shock or impact. In case of damage, sit down contact your dealer immediately to have the device checked for safety.
3. Impact or impact on the device may result in equipment failure as it is sensitive and highly accurate components are installed.

Precautions related to warranty

1. The manufacturer's warranty is 1 year from the date of purchase, the warranty on any material & workmanship 2 years.
2. The deadlines and conditions for a return policy of the merchant where you purchased our product apply.
3. The repair costs due to faults can be charged to the customer.

WARRANTY TERMS

The manufacturer guarantee is 1 year from the date of purchase.

For end customers in Germany and the EU, the warranty for material and workmanship defects is 2 years from the date of purchase.

CONDITIONS

The warranty period begins on the date of first sale of the product in new condition from a dealer to an end customer.

To assert warranty claims, the end customer must provide the original proof of purchase (i.e., the proof of purchase for the purchase of the product when new by an end user) and the warranty card.

EXCLUSIONS AND RESTRICTIONS

This Limited Warranty applies only to hardware components manufactured by or for Plabio that identify Plabio-initiated trademarks, trade names or the "PLABEAU" logo.

Plabio does not guarantee that the operation of the device is uninterrupted or error-free. Plabio is not liable for damages caused by failure to comply with the instructions for use were caused.

This warranty does not cover: (a) any damage resulting from the use of other products in combination with Plabio products were caused; (b) damages caused by accident, abuse, natural disasters or other external causes; (c) damage caused by the operation of the product outside the permissible or described by Plabio in the manual intended use were caused; (d) damages caused by the service (including repair service) which was performed by someone who is not a representative of Plabio; (e) for a device or part that has been modified without the written permission of Plabio to change its functionality or performance; (f) wearing parts, such as batteries, provided that no damage has occurred due to a defect in materials or workmanship; or (g) if a serial number has been removed or blurred by Plabio and the warranty card is lost.

OTHERS

✓ WEEE registration number PLABEAU: DE 79232390 (stiftung-ear)

This marking on the product, accessories, or literature indicates that the product and its electronic accessories (e.g., charger, USB cable) should not be disposed of with other household waste.

✓ WEEE registration number Lithium-Ion Battery: 12081500 (GRS Number, Umweltbundesamt)

This marking on the battery, manual or packaging indicates that the batteries in this product should not be disposed of with other household waste.

✓ Licensing of the packaging material used on: www.lizenzero.de

Information obligations according to §18 ElektroG (German electrical and electronic equipment law)

Disposal of electrical and electronic equipment

The crossed-out wheeled bin means that you have a legal obligation to transfer these devices to a separate collection from unsorted municipal waste. Disposal via the residual waste bin or the yellow bin is prohibited. If the products contain batteries or accumulators that are not firmly installed, they must be removed before disposal and disposed of separately as a battery.

Ways of returning old equipment (in Germany)

Owners of old devices can hand them in as part of the possibilities of returning or collecting old equipment set up and available by public waste disposal companies, in order to ensure the proper disposal of old equipment. Under the following link it is possible to have an online directory of collection and return points displayed:

<https://www.ear-system.de/ear-verzeichnis/sammel-und-ruecknahmestellen#no-back>

Privacy

We remind all end users of old devices that you are responsible for the deletion of personal data on the devices to be disposed.

WEEE registration number

Under the registration number DE 79232390 we are a part of the foundation elektro-altegeräte Register, Benno-Strauß-Str. 1, 90763 Fürth (Germany) registered as a manufacturer of electrical and / or electronic equipment.

Made in Korea

*** Read and save all instructions before using the product.**

Please direct any questions to the Plabio Europe Customer Hotline: **+49 6172 9258 127**
or by e-mail: **kontakt@plabio-europe.eu** or the local dealer.

Contact customer center for any inconvenience.

Please have your device model number, device status and contact number ready before contacting the customer center.

Visit the PLABEAU website or search the Internet for more comprehensive information.

www.plabio-europe.eu

Plasma changes your life

PLABIO

Manufactured by Plabio Co., Ltd. 10F, 316, Yeongdong-daero, Gangnam-gu, Seoul, Korea
Imported by Plabio Europe UG Hessenring 121, 61348 Bad Homburg, Germany

PLA BEAU
Plasma beauty solutions

G4+

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Deutsch 02
English 17